



TABLE OF CONTENTS

► [Home.md](#)

Home.md

Setup

Hardware

Physik Instrumente P-840.20 needs to be connected to the computer via USB.


Software

Install via PiezoaActuatorSetup.exe. Matlab Runtime 2024 needs to be installed. Physik Instrumente Software Suite could also be required.

How to:

Move

Connect the P-840.20 through the menubar and run the autozero function after connection has been established. Change the position or the voltage value to move the actuator accordingly. The plot documents the movement of the acuator.

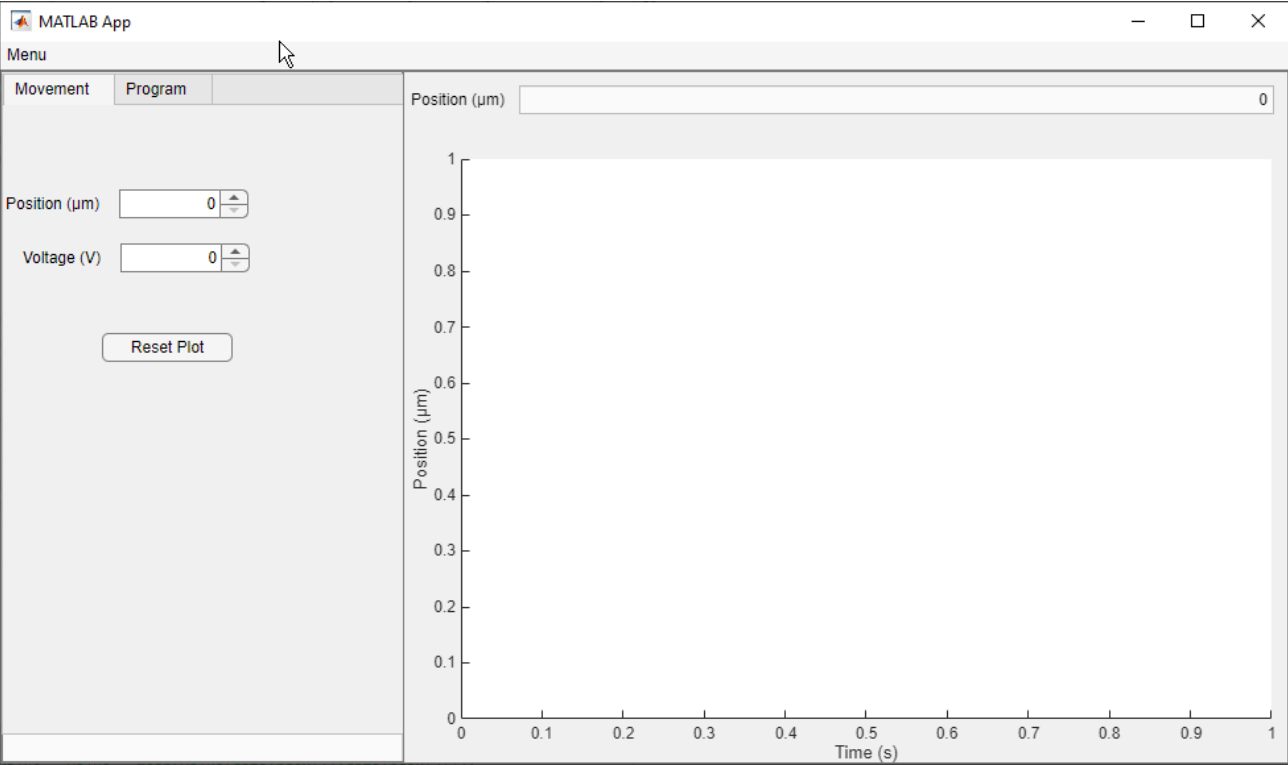

Program

Connect the P-840.20 through the menubar and run the autozero function after connection has been established. Change the various parameters to create the wanted movement protocol. The plot documents the movement of the acuator.

Menu

Movement

Program

Position relaxed (μm) Position extended (μm) Time relaxed (s) Time extended (s) Steps Speed ($\mu\text{m} / \text{ms}$) Time (ms)

Start

Stop

Ran AutoZero procedure

